



THUASNE



UniReliever™

Offloading single-upright
knee brace

UniReliever™

Universal size



TM5 Hinge ▶

(Lateral view)

Lightweight and easy to wear brace for knee osteoarthritis.

This universal OA bracing solution offers a targeted relief and personalized fit, with TM5 hinge technology.


Featuring correction targets and individualized force offloading, the cross-strap system with dial technology makes it easy for patients to control their correction for pain relief.

INDICATED FOR:

- Symptomatic unicompartmental femorotibial osteoarthritis (moderate to severe).
- Knee off-loading for post-traumatic, post-operative or degenerative conditions.
- Joint instability/laxity.
- Alternative to osteotomy or leg malalignment surgery.

ONLY 2 VERSIONS:

Left Medial / Right Lateral (U3090129900351)
or
Right Medial / Left Lateral (U3090129900352)

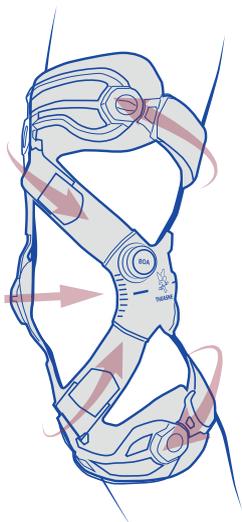

EFFICIENCY

- Dial fit system for precise offloading capacity.
- Combined cross strapping system with buttress designed to alleviate OA through even load distribution.
- TM5 hinge technology with range of motion control:
Extension limitations: 0°, 5°, 10°, 15°, 20° & 30°.
Flexion limitations (in option): 0°, 30°, 45°, 60°, 75°, 90° & 110°.

EASE OF USE



(Medial view)

Single dial for a more convenient correction experience combined with patella targets for easier alignment with knee center.



Quick release buckles with strap sequence labeling for easy donning.

COMFORT

- Padding with non-slip silicone webbing.
- Flexible open cuffs for a comfortable and breathable fit.
- Formable uprights and trimmable straps.

UniReliever™

Flexible cuffs adapt to most patients morphology

Consolidated strapping system with force buttress allows even disbursement of load at the point of application



Q1 Does the patient have chronic femoro-tibial OA⁽¹⁾ with a **dominant** compartment?



Q2 Does the patient fit within Thuasne standard size range?




Q3 Does the patient have a combined OA with knee instability?



Q4 What is the patient body mass index?



Q5 Does the patient have hand dexterity issues?

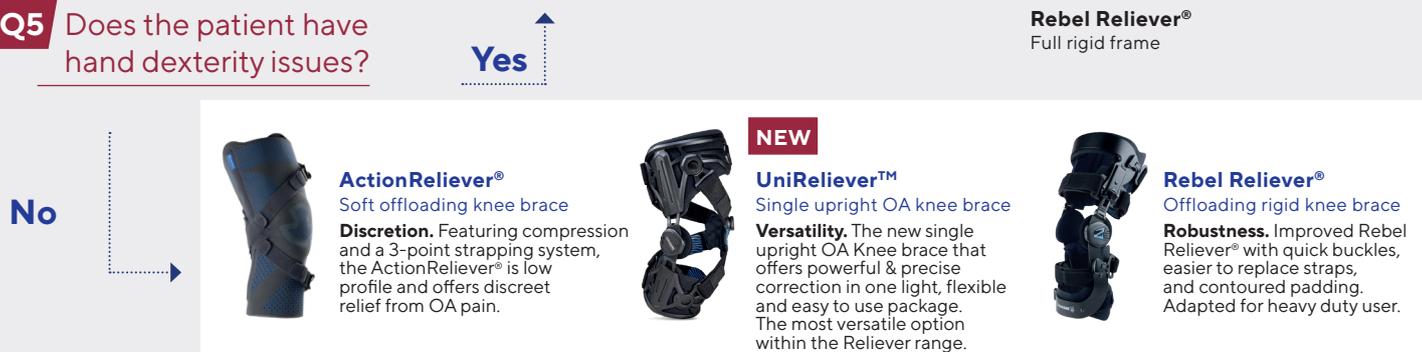

If the patient can be equipped with each of these products, please make sure to determine in **ACCORDANCE** with the patient, which product would suit the best his **EXPECTATIONS** and **NEEDS** (discretion, desired activities, ...).