

ABOUT COYOTE

At Coyote we are not only concerned with making quality products and providing excellent customer service, but also with creating a quality place to work and live. Our efforts toward a friendly, environmentally responsible workplace have earned us a Blue Ribbon Small Business Award from the U.S. Chamber of Commerce.

Made in the USA

More than 90 percent of our products are made in the USA. Most are made within a few miles of our home office in Boise, Idaho.

Need more product info?

Visit us at www.coyote.us for more information, videos, tips and instructions.







TABLE OF CONTENTS Warranty1

| Air-Lock2 | |
|-----------------------------|--|
| asy-Off Lock4 | |
| Grommet Lock6 | |
| Orop-In Locks7 | |
| Connectors8 | |
| Accessories10 | |
| anyard11 | |
| ocking Lanyard12 | |
| Proximal Lock14 | |
| Summit Lock16 | |
| Vave Suction Band17 | |
| Dynamic Strut AFO18 | |
| Small Pediatric AFO Strut20 | |
| Adhesive22 | |
| /alves24 | |
| ntegrator ML Gauge25 | |
| Coyote Composite26 | |
| ilent Air Compressor29 | |
| ools & Accessories30 | |
| ilent Air Compressor31 | |
| (leenaire32 | |
| ducationBack | |



2-Year Warranty

You want products to last and we want to learn how to make them last longer. This is why we offer a TWO YEAR warranty from date of purchase on all our products* against manufacturer's defects, excessive wear or breakage. You don't pay for shipping.**

- 1 Call and let us know if you need a part replaced. (208) 429-0026.
- We will need a PO to ship you a new part. We will then invoice you for the part and the shipping.
- 3 We will fax or email you a Return Authorization Form for you to complete.
- 4 Ship the damaged part and completed form to us. Charges bill to our UPS account.
- 5 After inspection and approval*** of returned part, we will issue you a credit. No charge for ground shipping.
- * A 90 day limited warranty applies to the Dynamic Strut AFO (CD206R, CD206M, CD206S) call for details.
- ** Parts must be used under reasonable and normal circumstances.
- *** All Kleenaire blowers are warranted for 1 year from date of original purchase.

 See warranty for details.

419 N. Curtis Rd. Boise, Idaho 83706 USA www.coyote.us (208) 429-0026



AIR-LOCK WASA

Coyote

A lightweight, dual-suspension pin system designed to combine suction suspension comfort with pin suspension security. The airtight seal prevents air from coming into the bottom of the socket, creating suction and reducing pistoning even without a sleeve! For patients who struggle with pin engagement, use the Deep Air-Lock with the 11 click pin.

- Suction allows for less traction on the limb for increased comfort and decreased liner wear
- Expulsion valve or suspension sleeve are not needed to maintain suction in most cases
- Water-resistant
- Do not use with elevated vacuum





CD103WH Lock wrench Makes maintenance easy.

www.coyote.us/airlock













CD103MD - Medium Deep Air-Lock



CD103SD - Small Deep Air-Lock



CD103AF **Alianable Connector** Offset up to .6"

ADDITIONS



One-Shot Connector Offset up to .6"



CD103SDI Single Direction Insert Offset up to .33"



CD103P8 10 mm, Two Part 8-click Pin Multi-Direction Insert 350 lbs capacity Offset up to .33" 1.3" long. Includes spacers







10 mm. 8-click Brass For active, heavier patients 350 lbs capacity 1.3" long. Includes spacers : 1.65" long. Includes spacers



CD103P11H 10 mm. 11-click Brass For active, heavier patients 350 lbs capacity



CD103PX 6 mm, Two Part 8-click Pin 1.3" long. Includes spacers



1.3" long. Includes spacers

ADDITIONS

EASY-OFF TH LOCK WSA

Our second generation lock has a more robust lever, tooth and larger O-ring. It is very similar to the Air-Lock in that it stays airtight to create suction in the socket, the Easy-Off Lock can be used with elevated vacuum. Its unique lever design also makes doffing easy, especially for amputees with hand strength issues, arthritis or trouble reaching the release button.

- Easy donning and doffing with lever release
- Self-lubricating materials reduce wear
- (>) Water-resistant
- (2) Can be used with elevated vacuum
- Airtight "suction" or vacuum suspension with the safety and security of a pin









Patent No. 6334876 Other patents pending



www.covote.us/easvoff



CD117D | Deep Easy-Off Lock

Easy removal without holding release button. Push or pull to release.



Pins are now sold separately and come in a two pack.



CD103AF **Alignable Connector** Offset up to 0.6"



CD111 One-Shot Connector Offset up to 0.6"



CD103SDI **Single-Direction Insert** Offset up to 0.33"



CD103MDI Multi-Direction Inser Offset up to 0.33"



CD103P8 10 mm, Two Part 8-click Pin 350 lbs capacity 1.3" long. Includes spacers



CD103P11 10 mm, Two Part 11-click Pin 350 lbs capacity



CD103P8H 10 mm, 8-click Brass 350 lbs capacity 1.65" long. Includes spacers: 1.3" long. Includes spacers: 1.65" long. Includes spacers



CD103P11H 10 mm, 11-click Brass For active, heavier patients: For active, heavier patients : 350 lbs capacity



CD103PX 6 mm, Two Part 8-click Pin 1.3" long. Includes spacers



1/4x20, Two Part 8-click Pin

1.3" long. Includes spacers

Coyote Catalog | 5

GROMMET LOCK WSA

Designed with pediatric patients in mind, the Grommet Lock is simple, lightweight and durable for use with upper and lower extremity. Just like our adult locks, it has multiple connector options for fabrication of a lower extremity prosthetic.

Grommet Lock





- Pediatric only
- (2) Lightweight, durable plastic
- (V) Water-resistant
- Ocan be modified for minimum clearance
- Can be used in upper and lower extremity prosthetics
- Three thread options for any liner type (10mm, 6mm, $\frac{1}{4}$ x 20)
- Weight limit of 100 lbs

Pins are now sold separately and come in a two pack.



CD104P10 10mm Pin



CD104P6 6mm Pin



CD104PO $\frac{1}{4}$ x 20 Pin



CD104P6H | 6mm Pin



www.coyote.us/grommet

CD104POH $\frac{1}{4}$ x 20 Pin

DROP-IN LOCKS WSA Coyote



CD122A | Air-Lock Drop-In

Build Height: 2.010"

First locks designed for 3D printed sockets and traditional fabrication.

- (Vater-Resistant
- Airtight suction dual suspension security of a pin
- Expulsion valve or suspension sleeve not needed in most cases
- (2) Great for 3D printed sockets For 3D printing contact Coyote® for your STL and OBJ Files.

CD122E | Easy-Off Drop-In

2nd Generation - Build Height: 2.175"



Easy donning and doffing with a more robust lever release and tooth

CD122AFD | Drop-In Air-Lock

Fabrication Dummy

CD122AFD - Fabrication dummy sold separately.



www.covote.us/dropin

CD122EFD | Drop-In Easy-Off Fabrication Dummy

2nd Generation

CD122EFD - Fabrication dummy sold separately.

Must use 2nd Generation with the current lock design



Coyote Catalog | 7

CONNECTORS WSA





CD103AF | Alignable Connector

Up to 0.6" of offset, 35a, Ideal for use with Air-Lock and Easy Off Locks, Can be glued to the lock in any offset position for fabrication into the socket, or keys directly onto locks when used with inserts. Eliminates need for offset components. For laminations or thermoforming, 350 lb weight limit,





CD103MDI | Multi-Direction Insert

Up to 0.33" offset in two directions, Installed into connector to allow connector to be keyed onto bottom of lock in centered position or offset in two directions. For laminations or thermoformina. Use with Alianable Connector or One-Shot Connector





CD103SDI | Single-Direction Insert

Up to 0.33" of offset in one direction. Installed into connector to allow connector to be keyed onto bottom of lock in centered position or offset in one direction. For laminations or thermoforming, Use with Alianable Connector or One-Shot Connector





CD111 | One-Shot Connector

Up to 0.6" offset, 44q. Connector with tie off ring for easier lay-ups. Lower profile than Alianable Connector. Can be alued to lock in desired location or keyed onto lock using inserts. For laminating only, 350 lb weight limit.





CD119SC | Test Socket Connector

Economical, cost effective, easy way to alue on connector designed specifically for Vivac check sockets. Connects to standard four hole componentry with our new high strength bridging design. Bonds well with Coyote Quik Adhesive. 11/16"height, 2 3/4" diameter, reusable.





CD106 | Alignment Coupler

Offset for ML and AP of 0.64". Easily installs on four hole connector. Easy to adjust with 3mm wrench. For fitting purposes only. 265 lb weight limit.





CD103PFF | Pediatric Fast Four Connector

Pediatric four-hole connector that keys directly onto the center of the lock for easy and fast fabrication. Used with Grommet Lock, Small Air-Lock or Medium Air-Lock. For laminations or thermoforming. 100 lb weight limit.



CD103PAF | Pediatric Alignable Connector

Pediatric four-hole connector that keys directly onto the center of the lock for easy and fast fabrication. Used with Grommet Lock, Small Air-Lock or Medium Air-Lock, Laminating or thermoforming. This model can be used to offset the connector from the bottom of the pin in virtually any location by up to 0.25".



CD115CF5 | 5 Degree AK Connector

Capable of up to 1.5" of Offset with 5 degrees of flexion angle. Connectors primarily use is with Lanyard Pin (CD118P) lanyard locks (CD103L, CD117L) can also be used with other Coyote® locks or the Integrator. Eliminates the need for offset components. For lamination or thermoforming, 265 lb, weight limit.

ACCESSORIES WSA





CD103FL | Fitting Lock

Air-Lock housing with reinforced distal end for accurate evaluation of lock pin spacing to improve seating of pin in lock. Ideal for reducing play in lock. Contact Coyote Design for long pin (Deep Air-Lock) fitting locks. For fitting only.

CD103RBC Release Button Cover

Added to end of release button to make it easier for patients to release the lock. Ideal for geriatric patients or others with lower dexterity or limited hand strenath.



CD103GPN | Guide Pin

Used to make standard locks deeper. Prevents long pins from bottoming out in standard locks.



CD103ELP **Exo Lamination Plug**

Used when fabricating the Air-Lock in Exo sockets.

CD103PC | Pin Cap

The Pin Cap is used to reseal socket to keep it airtight when deepening socket when using a long pin with a standard lock.



LANYARD





This is a simple, lightweight, ultra low build height lanyard lock that helps pistoning. Our new design for the metal puck has the connector posts built in for easier layups. The lanyard is segmented making easy length adjustments. They work well for simple visual donning and doffing and have K1 to K4 activity levels.



www.coyote.us/Lanyard-puck

Lanvard Aluminum

Dummy with Posts

CD121PD

CD120 | Lanvard Plastic Puck





CD121P | Lanvard Metal Puck with Posts

Integrated 4 hole connector







CD120L | Lanyard & low profile screw kit

24 inch strapping, (puck not included)

Lanyard Metal Puck

Integrated 4 hole connector





CD121D |

Aluminum Dummy

Covote Cataloa | 11


Lanyard

Covote Cataloa | 10

LOCKING LANYARD





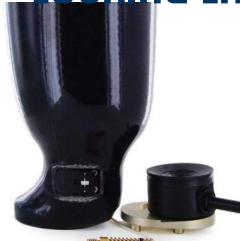

- The lanyard pin helps elongate the limb when donning.
- The pin gives audible confidence of engagement.
- The o-ring at the end of the pin creates an airtight socket helping to reduce air noise.

Fabrication with the 5 Degree AK Connector eliminates the need for offset components.

Same simple fabrication as the Deep Air-Lock and Deep Easy-Off Lock pin systems.

CD117L | Easy-Off Locking Lanyard

2nd Generation, Robust Lever Design / Push Button, lanyard lock

- Easy donning and doffing with lever release
- (>) Water-resistant
- Airtight "suction" with the safety and security of a pin





CD103L | Air-Lock Locking Lanyard

Push Button Design, lanyard lock

- Suction allows for less traction on the limb for increased comfort and decreased liner wear
- Expulsion valve or suspension sleeve are not needed to maintain suction in most cases
- () Water-resistant



CD115CF5 **5 Degree AK Connector**

Capable of up to 1.5" of Offset with 5 degrees of flexion angle.



CD118P **Locking Lanyard Pin**

350 lbs capacity, 1.65" long. Includes cord and spacers.



CD118PD **Lanyard Pinch Disk** Holds Lanyard cord to side of socket.

One pin is included with the locking lanyard locks

DDITIONS

PROXIMAL LOCK WSA

Simply easier. Designed to control rotation and be more user-friendly, the Proximal Lock attaches over a larger area, which distributes weight more evenly. The adhesive included with each lock is flexible even after it sets, making donning simpler.



The Proximal Lock now includes a reinforcement strip for extra reinforcement, less friction and repair.



and liner attachment.



CD116SA | Silicone Adhesive for Proximal Lock

Additional adhesive if needed. A tube is included with each lock.



- > ZERO CLEARANCE lock
- (Controls rotation
- (2) Use with your choice of distal adapters
- (V) Ideal for long limbs with minimal clearance
- \(\rightarrow\) Larger attachment distributes more tension
- Flexible attachment rolls on easier
- O Can be used with seal band liners
- Adhesive included with each lock
- (2) Lock can be retrofit
- Bills like a distal pin lock system (L5671)



www.coyote.us/proximal

CD116 | Proximal Lock









fabrication dummy, glue form, and adhesive.

Patent No. 6666894, 7431738, 7077868 Other patents pending.



SUMMIT LOCK WSA

Specifically for long limbs, the Summit Lock is designed to control rotation and requires zero clearance! There's no need to purchase a locking liner and you can use any socket adapter.

CD101 can be retrofit.



Base Lenath 1.75"







www.coyote.us/summit

- Liner attachment shape makes the liner simple to roll on
- Easy to cover cosmetically
- Ergonomic tab feeds easily into lock

CD102 | LowPro Summit Lock AK BK AE BE PED Patent No. 6666894

Base Lenath 2.5"



www.coyote.us/lowpro

Liner attachment is applied with "Coyote Quik Glue" or a 5 minute epoxy. Do not use silicone based adhesive.





CD402 | 20-Pack Extra Mix Tips For 50cc cartridae

WAVE SUCTION BAND

Wave Suction Band

The Coyote Wave Suction Band comes in multiple sizes.

- The Suction Band helps create dual suspension in combination with our Proximal, Summit and Low-Pro Locks.
- The airtight suspension produces a more intimate fit with added security.

The Hands-Free Valve or Push Button Valve are recommended to relieve suction when donning and doffing.

| Wave Suc | ave Suction Band sizes | | | |
|-------------|---|---|--|---|
| Part Number | Size | Circumference | Width | Thickness |
| CD130BS | Small | 16cm - 25cm | 7.62cm | 4mm |
| CD130BM | Medium | 25cm - 37cm | 7.62cm | 4mm |
| CD130BL | Large | 37cm - 49cm | 7.62cm | 4mm |
| CD130BX | Extra Large | 49cm - 64cm | 7.62cm | 4mm |
| CD130BXX | Extra Extra Large | 64cm - 77cm | 7.62cm | 4mm |
| | Part Number CD130BS CD130BM CD130BL CD130BX | Part Number Size CD130BS Small CD130BM Medium CD130BL Large CD130BX Extra Large | CD130BS Small 16cm - 25cm CD130BM Medium 25cm - 37cm CD130BL Large 37cm - 49cm CD130BX Extra Large 49cm - 64cm | Part Number Size Circumference Width CD130BS Small 16cm - 25cm 7.62cm CD130BM Medium 25cm - 37cm 7.62cm CD130BL Large 37cm - 49cm 7.62cm CD130BX Extra Large 49cm - 64cm 7.62cm |

Silicone adhesive is included.

The sizes are guidelines and can vary depending on the patient. Smaller sizes are more snua.





www.covote.us/wave-band

DYNAMIC STRUT AFO WSA

Our New Dynamic Strut Design combines a unique combination of strength and flex that other posterior dynamic struts do not offer. The AFO has the fit and comfort of a custom AFO with a posterior mounted dynamic strut for optimal ambulation. The strut is mounted to a custom thermoformed AFO and designed for average daily use; walking, working, hiking, biking, golf.

Has not been tested with running.

Designed to Flex

- Moves with the Patient
- Helps create a Natural Gate
- Works great with a thermoformed orthosis
- Offloads weight from the foot and ankle
- Provides Energy Response

Testimonials and demonstration videos can be found at www.coyote.us



www.coyote.us/dynamicstrut



AFO Strut Choices

CD206R (Rigid Flex)

Patient Weight Range: 175-265 lbs Patient Activity Level: Medium Activity

CD206M (Medium Flex)

Patient Weight Range: 125-200 lbs

Patient Activity Level: Light to Medium Activity

CD206S (Soft Flex)

Patient Weight Range: 125-200 lbs

Patient Activity Level: Light to Medium Activi-

CD206K (Strut Kit)

Kit includes Soft, Medium, and Rigid Flex

Struts and one set of parts.

Parts Included



A. Strut (1) (10.5 x 1.5)

- B. Screws (4)
- (2) 6x18 for Top Mounting Plate (2) - 6x20 for Bottom Mounting Plate
- C. Black Screws 6x14 (4)
- D. Clamp Washer (4)

Rigid Flex - CD206RS Medium Flex - CD206MS Soft Flex - CD206SS

Mounting Plates Kit - CD206MP

E. AFO Mounting Plate (2)

- (1) Upper Mounting Plate Short Posts
- (1) Lower Mounting Plate Tall Posts
- F. AFO Fab Tooling Pieces (2)

G. AFO Plate Foam (2)



Coyote

Counter

SMALL PEDIATRIC AFO STRUT





CD205PS | Small Pediatric AFO Strut

Patient Weight Range: up to 90 lbs. Patient Activity Level: All Activity Levels

A more streamlined AFO that vents well. Lighter weight than a solid AFO.

Practitioner Points

- Made from a composite plastic
- More torsion (Axial Rotation)
- Helps with hemiplegia, Cerebral Palsy, Spina Bifidg when there is correction needed in the foot.
- Helps with Tri-planar control along with control in the Sagittal and Coronal planes.
- Dynamic return with Plantar flexion Resist/ Dorsiflexion Resist



Technical Points

- Simple design
- Adds strength, reinforcement and support to Thermoplastic AFOs.
- All activity levels for children up to 90 lbs.

(Various Locking concave design for strength and thermoplastic vacuum drawing compatibility.







ADHESIVE WSA





Ideal Multi-Purpose Adhesive

Great for any situation when you need something glued quickly and easily. Both Adhesives ship nonhazardous and are safe with no odor.

Very quick set with no sag



CD4220 **Quik Glue White**

30 Second, 220 cc cartridge. Dispenser required; includes 1 mix tip.

30 Second set up time 2 minute full cure time

CD4150 **Ouik Glue White**

30 Second, 50cc cartridge. Dispenser or plunger required; includes 1 mix

NEW ADHESIVE



(2) 60 second with no sag

- () Is solid when cured
- (2) Bonds well to fabric, plastic and metal
- () Easy to smooth

CD4320 **Quik Glue Tan**

60 Second, 220 cc cartridge. Dispenser required: includes 1 mix tip.

60 Second set up time 3 minute full cure time

CD4350 Ouick Glue Tan

60 Second, 50cc cartridge. Dispenser or plunger required; includes 1 mix tip.

Stays Flexible when cured



- () Longer set time with no sag
- Stays flexible when cured
- (2) Bonds well to fabric, plastic and metal
- () Easy to smooth

Coyote **QUIK GLUE** Coyote SET UP TIME **OUIK GLUE** BLACK - CD4220S 7 minute full cure time - 220cc BLACK-CD4250 All purpose urethone adhesive designed to set up fast, hold shape with no sag

CD4220S Ouik Glue Black

90 Second, 220 cc cartridge. Dispenser required; includes 1 mix tip.

90 Second set up time 7 minute full cure time

CD4250 **Quick Glue Black**

90 Second, 50cc cartridge. Dispenser or plunger required; includes 1 mix tip.

^{*}Recommend purchasing additional tips - Tips must be replaced prior to each use.



Coyote

QUIK GLUE

WHITE - CD4220

2 minute full cure time : 220cc

to set up fast, hold shape with no sag-

VALVES WSA



Coyote Design's auto expulsion valves are easy to install in new sockets with fabrication tooling or retrofit. The built-in auto expulsion quietly expels air all day with no fuss and features a higher "crack" pressure than most valves, which helps to create a higher level of hypobaric suspension.

CD105 | CQL Auto Expulsion Valve

Height: 11mm. Patent pending.

- Easy to clean and retrofit
- Quiet expulsion
- () Low profile
- Includes a filter to block debris



CD105HF | HandsFree Valve

Height: 8.5mm. Patent pending.

- Works with seal band liners
- Manual twist open airflow
- Flexible inner style valve



CD103WH Lock wrench Makes maintenance easy.





www.coyote.us/valves

THE INTEGRATOR



The Integrator Socket Adapter, when used with one of our connectors, can be built into alignment in one step with a clean inside socket finish. Ideal for AK or BK sockets.



Integrator™ (2 3/8" diameter) AK BK





www.coyote.us/integrator

Standard size, fits most lower extremity limbs. Integrator glues to mold. Alignable Connector with or without inserts is keyed onto Integrator and both parts are fabricated into socket. 265 lb weight limit. Thermoforming or laminating.



Small Integrator™ (2"diameter) AK BK





Smaller size for smaller limb shapes. 265 lb weight limit. Thermoforming or laminating.



Put it to the Test

The Proximal lock, LowPro and Summit Locks were designed for Coyote® CEO and Co-Founder


Matt won 5 ITU Triathlon World Championships in 2006, 2007, 2008, 2009, 2010 while wearing the Proximal Lock.

COYOTE COMPOSITE

A great alternative to itchy, expensive carbon fiber braid! Coyote Composite, is made from basalt filament, a safe, tough material for composite sockets and AFOs. Also relatively easy to cut and finish compared to other materials.

- Extremely durable
- Resistant to cracking
- (Less itch than carbon
- Non-carcinogenic and non-toxic

- Strong and lightweight
- Finishes smooth, with superior resin saturation
- (2) Unique dampening characteristics
- (Less expensive

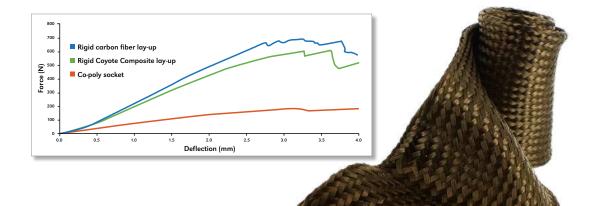



Resin Coyote Composite works with typical prosthetic resins and saturates much better than carbon.

Saturation and Laminating Because of its superior saturation it's important to prevent resin starvation. Lamination starvation will cause frayed edges.

Finishing Coyote Composite finishes easily. Use sharp tooling at slow speeds to finish. In some cases, you may need to wet-sand.

Relieving Coyote Composite is heat- and flameresistant and isn't damaged by heating. It allows for limited adjustment based on the type of resin used.

Safety Because of the inert nature of basalt, it is not a carcinogen. Also, Coyote Composite is between 9 and 16 micrometers in diameter, which is too large to be an inhalant risk. Basalt also causes significantly less skin irritation.

Composite on outside of lay-up



Composite on inside of lav-up



www.coyote.us/coyotecomposite

COYOTE COMPOSITE



Braid Our proprietary braid was designed to be used in equivalent amounts to carbon. Coyote Composite is more flexible than carbon (we consider this to be an advantage), but if you need to match the rigidity of a carbon socket, use additional layers of stockinette and/or carbon tape.



Fabric Coyote Composite fabric sheeting can be used for reinforcement layers of heavy duty sockets and/or AFOs. The material is thin and lightweight which allows for greater flexibility in adding the ideal amount of reinforcement

| abric | | Туре | TPI* | Length/ (ft) |
|-------|----|------|------|-----------------|
| ŭ | CD | F | 350 | 3 |

Fabric available in

.17 mm Thickness Has a Bi-directional weave

Ordering options

| | Туре | Size (in) | Length/ (ft) |
|----|------|-----------|-----------------|
| CD | В | 3 | 10 |
| CD | В | 3 | 20 |
| CD | В | 3 | 82 |
| CD | В | 4 | 10 |
| CD | В | 4 | 20 |
| CD | В | 4 | 82 |
| CD | В | 5 | 10 |
| CD | В | 5 | 20 |
| CD | В | 5 | 82 |
| CD | В | 6 | 10 |
| CD | В | 6 | 20 |
| CD | В | 6 | 82 |
| CD | В | 7 | 10 |
| CD | В | 7 | 20 |
| CD | В | 7 | 82 |
| CD | В | 8 | 10 |
| CD | В | 8 | 20 |
| CD | В | 8 | 82 |
| CD | В | 10 | 10 |
| CD | В | 10 | 20 |
| CD | В | 10 | 82 |

SILENT AIR COMPRESSOR



These silent air compressors are dependable, vibration free, and compact that have noise levels as low as 30 dB/A.



Each high performance compressor is built with quality in mind, and comes fully equipped with powder coated air tank, pressure switch, 5 micron line filter, pressure regulator, and gauges for easy installation and trouble free operation. When operating under full load this line of exceptionally quiet air compressors offers decibel levels as low as 30 dB/A.

Visit our website for more information of the different silent air compressor packages.



www.coyote.us/compressor

TOOLS & ACCESSORIES WEEK







Extractor Socket Removal Tool

The Extractor quickly and easily removes most molds from the laminated or thermoformed socket without damaging the mold or socket. Includes adapter for sockets with fourhole connectors

Modification Handle

Allows for modifications while keeping fingers and knuckles clear of abuse. Easy to clean and replace blades, and reduces blade breakage. Works with 10" Stanley Surform Blades, flat- or half-round

Plaster Spreader CD307

Used for build-ups during modifications. The tool is easy to clean and won't stick to plaster.



www.coyote.us/tools



12 oz Mixina Bowl


Made of non-stick, durable material, its easy to grip and clean. Wide base won't tip over.



Small Casting Dots (200 pk) Large Casting Dots (50 pk)

Easily mark problem areas during the casting process. The adhesive back will keep them in place during casting. 0.375" or 0.87" diameter.



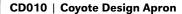

Take care of your clothes in the lab with this inexpensive Coyote Design Smock. One size fits all.

Casting Bag

The Casting Bag, as seen in Q School workshops, is ideal for casting over to keep the liner clean. Sold in packs of 100.

CD317 | ML Gauge

This ML Gauge was designed specifically for measuring prosthetic and orthotic patients. Priced right and built to last.









New Price - New Design

KLEENAIRE | Dual Filtered Air Purification System


Patent US D536, 432S

- Reduces spread of viruses
- Lowers CO₂ levels
- Easy to install, maintain, and cost effective
- Ocontrols dust, removes odors, contaminates, and even viruses
- Higher air quality and healthier working environment
- Designed for minimum noise and high air flow

The affordable alternative to purify the air for residential, comercial, or manufacturing





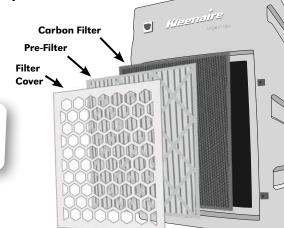

Coyote

Buy your **filter subscription** at www.coyote.us/kleenaire/

See Manufacturer's Capacity Index chart for information on contaminates removed by the charcoal filter.

For more information go to www.coyote.us







Coyote has education opportunities everywhere

www.coyote.us ● YouTube ● opqschool.com ● Webinars ● Conferences

Tell your colleagues about our continuing education credits.



